


# MetPrep Preparation Procedure – No 92

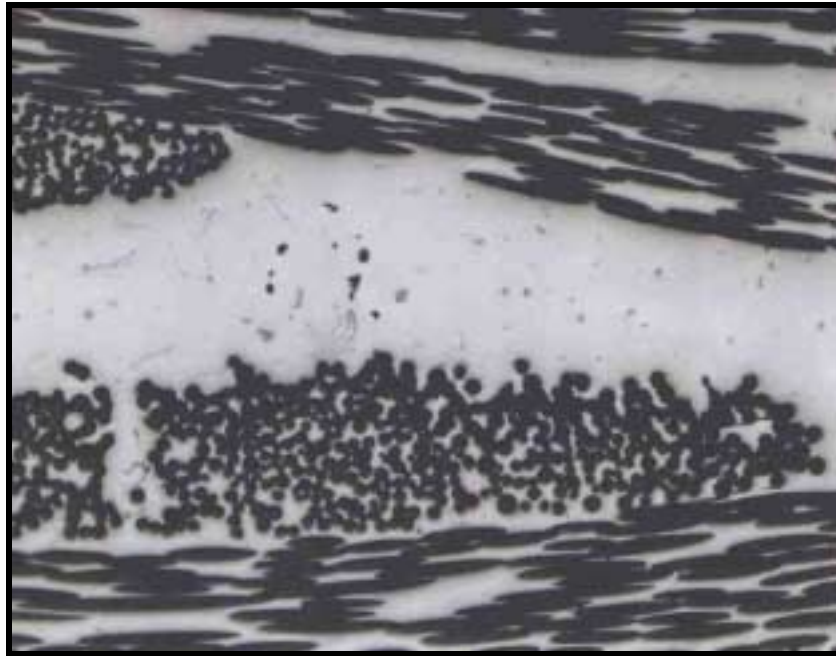

## Silicon Carbon fibres in Aluminium Alloy

	Surface	Abrasive	Pressure		Speed – Dir	Time
			Psi	N		
Primary Grinding Stage	Paper	P240g SiC	5	25	150 – Comp	Until Planar

	Surface	Abrasive	Pressure		Speed – Dir	Time
			Psi	N		
Additional Grinding Stages	Planocloth H	9 µm (WB) Diamond	5	25	250 – Comp	5 mins
	Planocloth	3 µm (WB) Diamond	5	25	250 – Comp	3 mins
	Durasilk	0.06 µm Silco	5	25	100 – Comp	3 mins

	Surface	Abrasive	Pressure		Speed – Dir	Time
			Psi	N		
Polishing Stage	Multicloth	0.06 µm Silco	3	15	80 – Comp	2 mins